

## **VERSATIP**

The versatile  
solution  
for Endo-Urology



VERSATIP



## The next generation instruments in Endo-Urology

- ▼ Detachable – one handle for all inserts
- ▼ Color code – the easy identification of the tip
- ▼ Two handles – self holding and standard
- ▼ Easy connection – the reliable ball tip technique
- ▼ Light weight handle – the perfect fit for your fingers
- ▼ Make cleaning easy – flushing port for the complete instrument
- ▼ 360° rotation - just with your finger tip
- ▼ All parts can be ordered separately



Watch the  
product video




360°



Multiple tips design



Color coded inserts

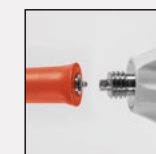

Unique, detachable  
system



Flushing port  
for easy cleaning




Two handles -  
multiple options

FORCEPS & PUNCHES

- Automatic self-holding
- Black color coded




351-100-550  
Self-holding handle for  
forceps and punches  
(recommended handle)

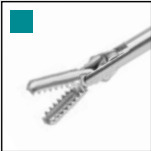

| Grasping Forceps with toothed cup,<br>double action, 600mm working length                      |                   | Insert only | Complete instrument |
|------------------------------------------------------------------------------------------------|-------------------|-------------|---------------------|
|  semi-rigid | 3 Charr., Ø 1.0mm | 253-950-503 | 253-950-003         |
|                                                                                                | 4 Charr., Ø 1.3mm | 253-950-504 | 253-950-004         |
|                                                                                                | 5 Charr., Ø 1.6mm | 253-950-505 | 253-950-005         |
|  flexible   | 3 Charr., Ø 1.0mm | 253-951-503 | 253-951-003         |
|                                                                                                | 4 Charr., Ø 1.3mm | 253-951-504 | 253-951-004         |
|                                                                                                | 5 Charr., Ø 1.6mm | 253-951-505 | 253-951-005         |



| Stone Grasping Forceps for large stones,<br>double action jaws, 600mm working length           |                   | Insert only | Complete instrument |
|------------------------------------------------------------------------------------------------|-------------------|-------------|---------------------|
|  semi-rigid | 3 Charr., Ø 1.0mm | 253-952-503 | 253-952-003         |
|                                                                                                | 4 Charr., Ø 1.3mm | 253-952-504 | 253-952-004         |
|                                                                                                | 5 Charr., Ø 1.6mm | 253-952-505 | 253-952-005         |
|  flexible   | 3 Charr., Ø 1.0mm | 253-953-503 | 253-953-003         |
|                                                                                                | 4 Charr., Ø 1.3mm | 253-953-504 | 253-953-004         |
|                                                                                                | 5 Charr., Ø 1.6mm | 253-953-505 | 253-953-005         |

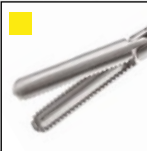

| HAMADA long Stone Grasping Forceps,<br>double action jaws, 600mm working length                  |                   | Insert only | Complete instrument |
|--------------------------------------------------------------------------------------------------|-------------------|-------------|---------------------|
|  semi-rigid | 3 Charr., Ø 1.0mm | 253-954-503 | 253-954-003         |
|                                                                                                  | 4 Charr., Ø 1.3mm | 253-954-504 | 253-954-004         |
|                                                                                                  | 5 Charr., Ø 1.6mm | 253-954-505 | 253-954-005         |
|  flexible   | 3 Charr., Ø 1.0mm | 253-955-503 | 253-955-003         |
|                                                                                                  | 4 Charr., Ø 1.3mm | 253-955-504 | 253-955-004         |
|                                                                                                  | 5 Charr., Ø 1.6mm | 253-955-505 | 253-955-005         |



| Biopsy Forceps, 600mm working length                                                             |                   | Insert only | Complete instrument |
|--------------------------------------------------------------------------------------------------|-------------------|-------------|---------------------|
|  semi-rigid | 3 Charr., Ø 1.0mm | 253-956-503 | 253-956-003         |
|                                                                                                  | 4 Charr., Ø 1.3mm | 253-956-504 | 253-956-004         |
|                                                                                                  | 5 Charr., Ø 1.6mm | 253-956-505 | 253-956-005         |
|  flexible   | 3 Charr., Ø 1.0mm | 253-957-503 | 253-957-003         |
|                                                                                                  | 4 Charr., Ø 1.3mm | 253-957-504 | 253-957-004         |
|                                                                                                  | 5 Charr., Ø 1.6mm | 253-957-505 | 253-957-005         |



SCISSORS

- Standard handle
- White color coded




351-100-500  
Standard handle  
for scissors  
(recommended handle)



| Scissors, single action,<br>600mm working length                                               |                   | Insert only | Complete instrument |
|------------------------------------------------------------------------------------------------|-------------------|-------------|---------------------|
|  semi-rigid | 3 Charr., Ø 1.0mm | 253-958-503 | 253-958-003         |
|                                                                                                | 4 Charr., Ø 1.3mm | 253-958-504 | 253-958-004         |
|                                                                                                | 5 Charr., Ø 1.6mm | 253-958-505 | 253-958-005         |
|  flexible   | 3 Charr., Ø 1.0mm | 253-959-503 | 253-959-003         |
|                                                                                                | 4 Charr., Ø 1.3mm | 253-959-504 | 253-959-004         |
|                                                                                                | 5 Charr., Ø 1.6mm | 253-959-505 | 253-959-005         |

ASSEMBLY

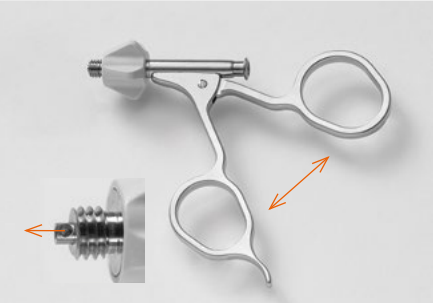

**STEP 1**

Open the handle.  
The ball socket pushes out.


Keep the handle open.



**STEP 3**

Push back the locking nut.


With both opened handle and jaws, connect the ball head and ball socket.

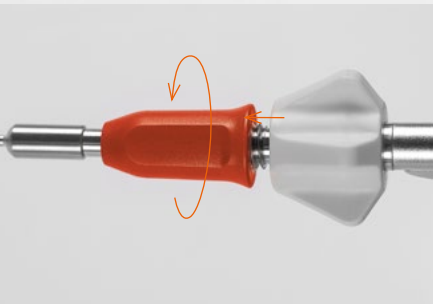

**STEP 5**

Close the handle. Shaft and insert will be retracted.  
Keep the handle closed.


Screw the locking nut onto the handle.

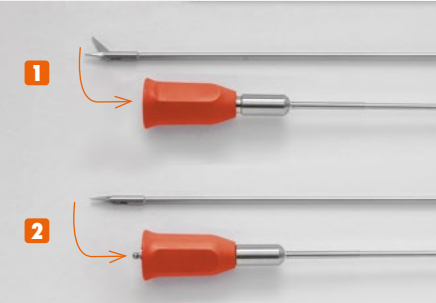

**STEP 2**

**1** Opened jaw

**2** Close the jaw - the ball head pushes out. Keep the jaw closed.

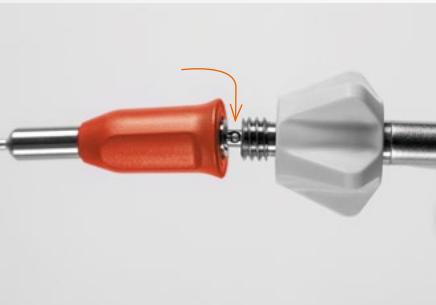

**STEP 4**

Insert the ball head into the socket.

Keep the handle opened.



**DONE**

Check the connection for being tight.

The knob rotates 360° and the handle as well as the jaw opens and closes.



RZ Medizintechnik GmbH

Unter Hasslen 20/22  
78532 Tuttlingen, Germany

+49 7462 9470-0

[sales@rz-medizintechnik.com](mailto:sales@rz-medizintechnik.com)  
[www.rz-medizintechnik.com](http://www.rz-medizintechnik.com)

